
PBG-866-IG-E
2-4/00hk
PK-IB-PIL-866-01-3

www.zepter.com

© COPYRIGHT HOME ART & SALES SERVICES AG - ZEPTER GROUP - Printed in Italy 2011 Code **PMD-000-00-ML** Artwork: BRID MARKETING - ITALY

 **Vital System**
by Zepter Group

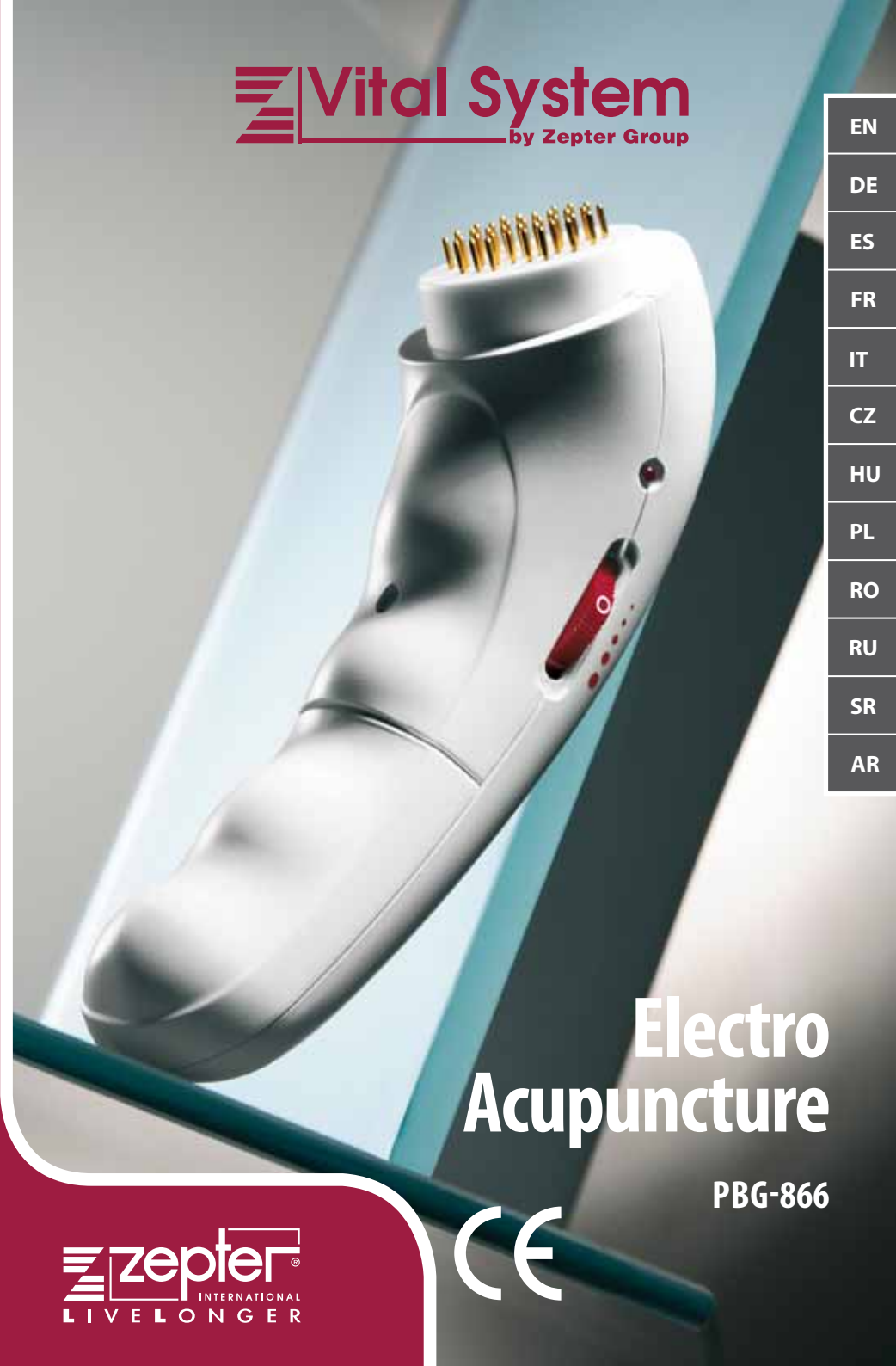EN
DE
ES
FR
IT
CZ
HU
PL
RO
RU
SR
AR

Electro Acupuncture

PBG-866

Index

LG-866

EN	Instruction Manual	2
DE	Betriebsanleitung	12
ES	Manual de instrucciones	22
FR	Manuel d'Instruction	32
IT	Manuale di istruzioni	42
CZ	Návod k obsluze	52
HU	Használati útmutató	62
PL	Instrukcja obsługi	72
RO	Manual de instrucțiuni	82
RU	Инструкция по эксплуатации	92
SR	Uputstvo za upotrebu	102
AR	Instruction Manual	112

Electro Acupuncture

Relieves Pain and Tension. Revitalizes and Invigorates the Body

FOR THE RELIEF OF SYMPTOMATIC PAIN ONLY. NEITHER A CURE NOR A MEDICAL DEVICE.

What is electro-acupuncture?

Welcome to an innovative approach to pain and tension management that couples an ancient method with space age technology.

ELECTRO-ACUPUNCTURE is a hand-held, very easy to use, do-it-yourself, battery-operated electronic device for self-relief through direct stimulation of pain, aches, tensions and stress-related insomnia. This stimulation can also serve to revitalize, invigorate and enhance the body's energy level.

ELECTRO-ACUPUNCTURE can go any where, with any one, at any time.

How was ELECTRO-ACUPUNCTURE Developed?

ELECTRO-ACUPUNCTURE was developed over a period of 20 years for use by the Soviet cosmonauts for simple, uncomplicated, immediate self-relief from pain, aches and tension. After many years of successful application in the Soviet space program, the device was further perfected and adapted for civilian use by a leading Russian scientist.

What is ELECTRO-ACUPUNCTURE Based On?

ELECTRO-ACUPUNCTURE is based on modern, space age technology coupled with the 4000 year old Chinese health care science of acupuncture.

Acupuncture has throughout the ages, been used, among other things, to relieve pain, aches and tension by direct stimulation of specific points in or near the nerves under the surface of the skin. In classic acupuncture the points follow a series of life-forcelines called 'meridians', which circle the body longitudinally. In this approach, the relevant point or points of a specific pain, ache or tension source must be accurately located by a trained acupuncturist, who then stimulates the point(s) by penetrating the surface of the skin with special needles.

Introduction

How is ELECTRO-ACUPUNCTURE Different?

ELECTRO-ACUPUNCTURE is a completely noninvasive, state-of-the-art method that replaces the classic acupuncture needle with electronic pulse stimulation. Stimulation is carried out with 36 goldplated, conductive pulse-pins. In addition, the device's unique BIO-FEEDBACK capability solves the problem of correct application by locating the exact point to be stimulated. The user simply applies the device to the overall area as indicated on the body maps (see pages 10-11). The device identifies the exact spot by detecting the low resistance points that are a critical mechanism in organic and natural pain relief. The pulse-pins then send an electronic signal to the identified spot on the body.

How is ELECTRO-ACUPUNCTURE Applied?

As indicated above, all the user needs to do is place the device on the area where pain, aches or tension are felt, or if preferred, on the classic acupuncture point (see drawings), turn it on, and let **ELECTRO-ACUPUNCTURE** go to work! The device locates the correct spot - modern or classic - to stimulate, while the user determines the intensity of the level of stimulation and the length of time applied. This self-relief application is totally natural, safe and with absolutely no sideeffects!

Which Pains and Aches Can Be Relieved with ELECTRO-ACUPUNCTURE?

ELECTRO-ACUPUNCTURE is effective in providing relief for a wide variety of pain conditions, including migraine and tension headaches, aching muscles and muscle strain, stress and stress-related insomnia, back, shoulder, elbow and knee pain, rheumatic, arthritic and joint degeneration related-pain, and more. It can also help revitalize and invigorate the body! While the symptomatic relief of some pains and aches can be achieved immediately, others respond to the cumulative effect of several applications.

The ELECTRO-ACUPUNCTURE device is not a cure and makes no claims to that effect.

WARNING: Application of **ELECTRO-ACUPUNCTURE** may be hazardous to the user under certain conditions, including, but not limited to those listed:

* **Do not apply the ELECTRO-ACUPUNCTURE** device on or near an open wound or wounds, lesions or eruptions resulting from skin disorders.

* **Do not use the ELECTRO-ACUPUNCTURE** device if you have a heart disease or any coronary problems, or have suffered thrombosis of any kind, or if you are pregnant.

In case of any doubt concerning the suitability of an **ELECTRO-ACUPUNCTURE** application, the user's health care service provider must be consulted.

Specific recommendations for the correct use of the **ELECTRO-ACUPUNCTURE** device are detailed on the following pages.

NOTE: The **ELECTRO-ACUPUNCTURE** device should be kept out of the reach of children.

HOW TO USE ELECTRO-ACUPUNCTURE?

Method 1: Application to General Area

With this method the user simply places the **ELECTROACUPUNCTURE** device on the general area of the pain, ache or tension, and begins application according to the following procedure:

1. Turn the device on with the on/off button on the side.
2. Temporarily select a pulse intensity level by turning the regulator knob to low, while holding the **ELECTRO-ACUPUNCTURE** device in your hand.
3. Place the pulsating **ELECTRO-ACUPUNCTURE** device on the area of your body where you feel pain, aches or tension (carefully avoiding the areas and conditions noted above in the **WARNING** section), and stimulate the selected area with rotating motions.
4. Now determine the correct pulse intensity for you by slowly turning the regulator knob, stopping when you feel a tickling sensation. This sensation, which should be comfortable for you to tolerate, indicates that you have arrived at the correct level of intensity for you. The intensity can be increased or decreased during an application in response to a change in sensation. It is recommended to use the high intensity level only with circular motions.
5. Stimulate the selected area with the device, preferably using rotating motions lightly over the skin, for 3 to 7 minutes (longer application time is possible if desired). As noted, it is not necessary to apply the **ELECTRO-ACUPUNCTURE** directly to the exact spot since the device will seek and find the correct point. Do not press the **ELECTRO-ACUPUNCTURE** deep into the skin as it may cause irritation.
6. If an **ELECTRO-ACUPUNCTURE** application is begun at the onset of a specific pain, ache or tension, relief is generally achieved following just one application. Repeat applications are recommended if the pain returns. In cases of longterm and/or chronic pains, better relief response is achieved on a cumulative basis, with a gradually increasing period of time between each application.

NOTE: A slight redness may appear during or after use but should disappear after a short period of time.

Method 2: Application to Classic Acupuncture Points


With this method the user places the **ELECTRO-ACUPUNCTURE** device on or near a specific acupuncture point, which, according to the ancient science of acupuncture described above, controls the well being of a certain part of the body. As can be seen in the drawings on the following pages, the classic acupuncture point may, in some cases, not even be near the part of the body that is in pain. If this method is chosen, the user should study the relevant drawings prior to selecting the classic acupuncture point to be used for an application. The procedure for application according to Method No.2 is the same as that described above for Method No.1. Note that with Method No.2 it is also not necessary to locate the precise classic acupuncture point. The user needs only to apply the device to the general area of the selected classic acupuncture point leaving the location of the exact spot to **ELECTROACUPUNCTURE** device's **BIO-FEEDBACK** capability.


**Relieves Pain and Tension
Revitalizes and Invigorates the Body**

Classic Acupuncture Points

FRONTAL-POINTS

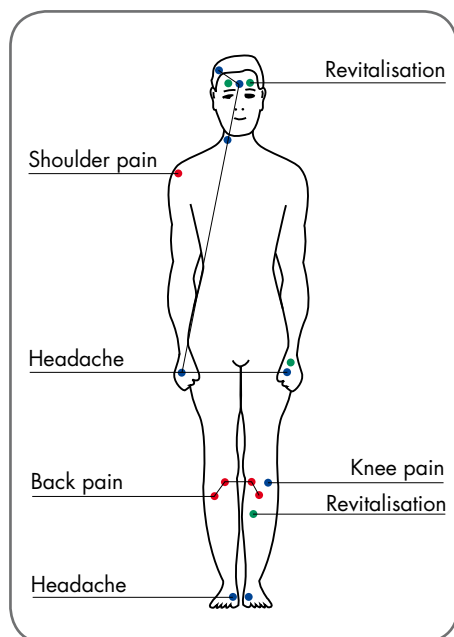BACKSIDE-POINTS

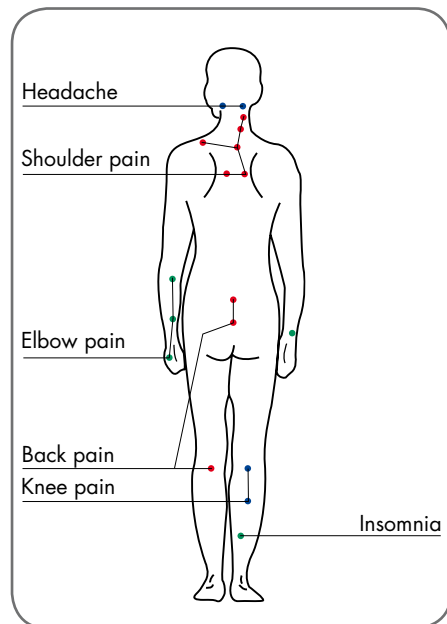Headaches - soothe your pain

Headaches have numerous causes, including migraine, congestion and sinus conditions, toothaches, over exertion, fatigue, tension, pressure, etc. Using a circular motion, stimulate the points shown in the drawings with the **ELECTRO-ACUPUNCTURE** device, for relief from headaches, carefully avoiding the areas and conditions noted in the WARNING section.

Fig.1
is most effective for
lower head pain.

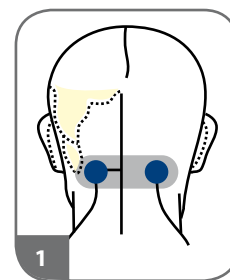
Fig.2
is most effective for
congestive conditions.






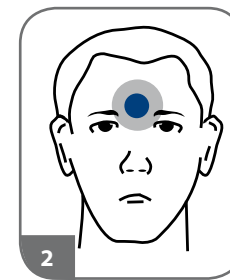
Fig.4
is effective for general
headaches.

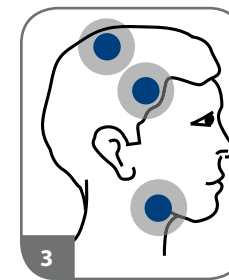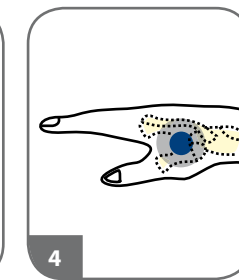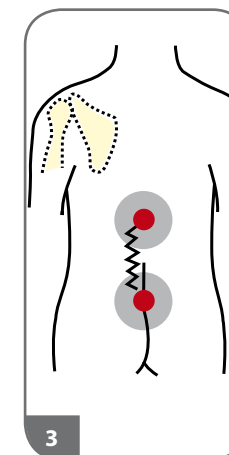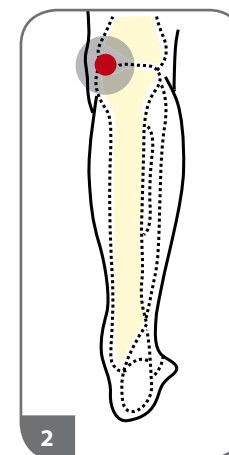Back-Pain - ease your pain

The most common back pain is in the lower back. Back pain can be caused by a sudden strenuous movement, bad posture, a pulled muscle, lifting incorrectly, menstrual pain, etc. Using a circular motion, stimulate the points shown in the drawings with the **ELECTRO-ACUPUNCTURE** device, for the relief of back pains, carefully avoiding the areas and conditions noted above in the WARNING section.

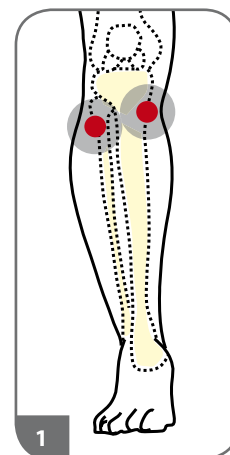Insomnia - Self Relief

Insomnia can be caused by tension or stress, change in sleeping environment, exposure to high altitude, or a sudden sensitivity to light noise, large meals, extensive exercise etc. Using a circular motion, stimulate the points shown in the drawings with the **ELECTRO-ACUPUNCTURE** device, for the relief of insomnia, carefully avoiding the areas and conditions noted above in the **WARNING** section.


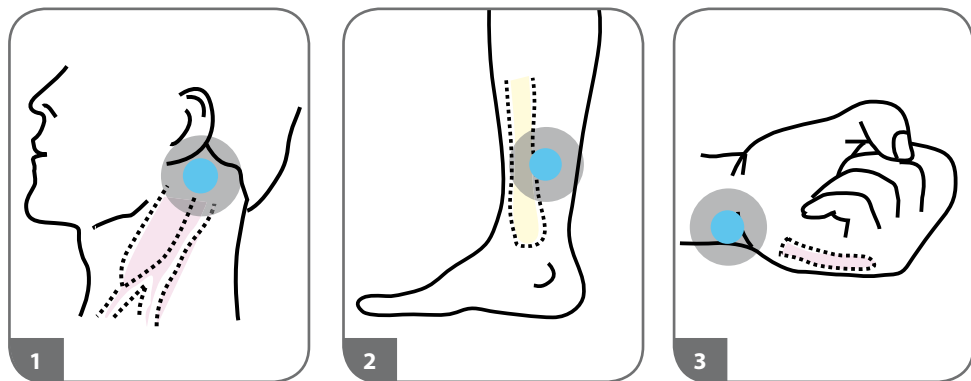
Fig.3
is the most effective for relief
of insomnia.

Shoulder-Pain - Relieve Your Pain

Shoulder pain can be caused by a pulled or strained muscle, over exertion, trauma or a direct blow to the shoulder, bursitis, arthritis, rheumatism, etc. Using a circular motion, stimulate the points shown in the drawings with the **ELECTRO-ACUPUNCTURE** device, for the relief of shoulder pain, carefully avoiding the areas and conditions noted above in the **WARNING** section.

Fig.2
is the most effective
for the relief of
shoulder pain.

Fig.4
is the secondary relief
points for shoulder
pain.

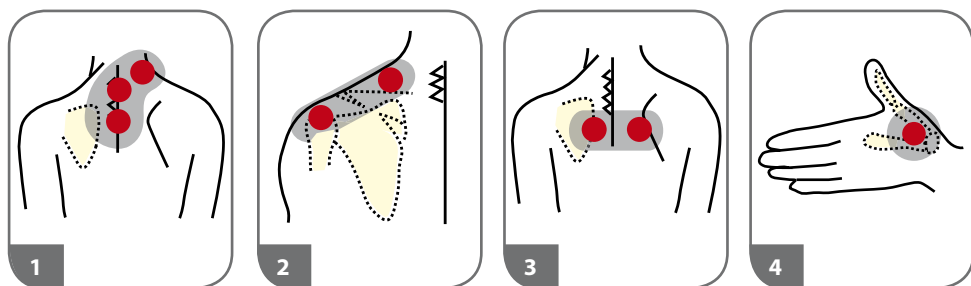Elbow-Pain - Calm Your Pain

Elbow pain can be caused by a pulled or strained muscle, over exertion, trauma or direct blow to the elbow, tendinitis, bursitis, arthritis, rheumatism, etc. Using a circular motion, stimulate the points shown in the drawings with the **ELECTRO-ACUPUNCTURE** device, for the relief of elbow pain, carefully avoiding the areas and conditions noted above in the **WARNING** section.

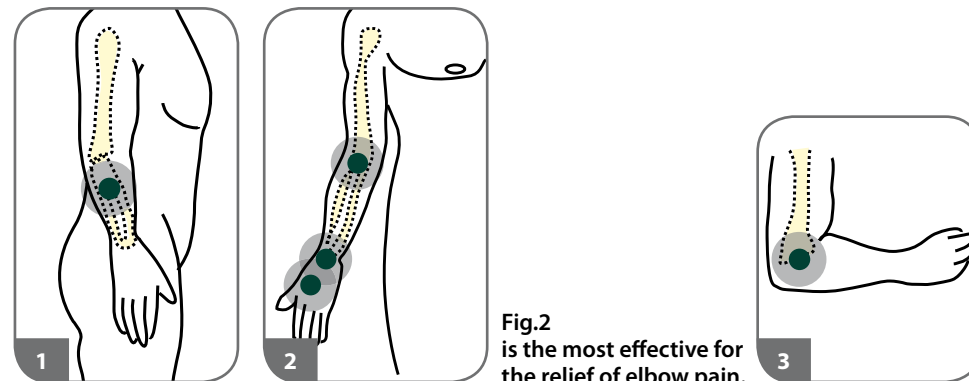
Fig.2
is the most effective for
the relief of elbow pain.

Knee-Pain - Alleviate Your Pain

Knee pain can be caused by a pulled or strained muscle, trauma or direct blow to the knee, arthritis, rheumatism, over exertion, etc. Using a circular motion, stimulate the points shown in the drawings below with the **ELECTRO-ACUPUNCTURE** device, for the relief of knee pain, carefully avoiding the areas and conditions noted above in the **WARNING** section.

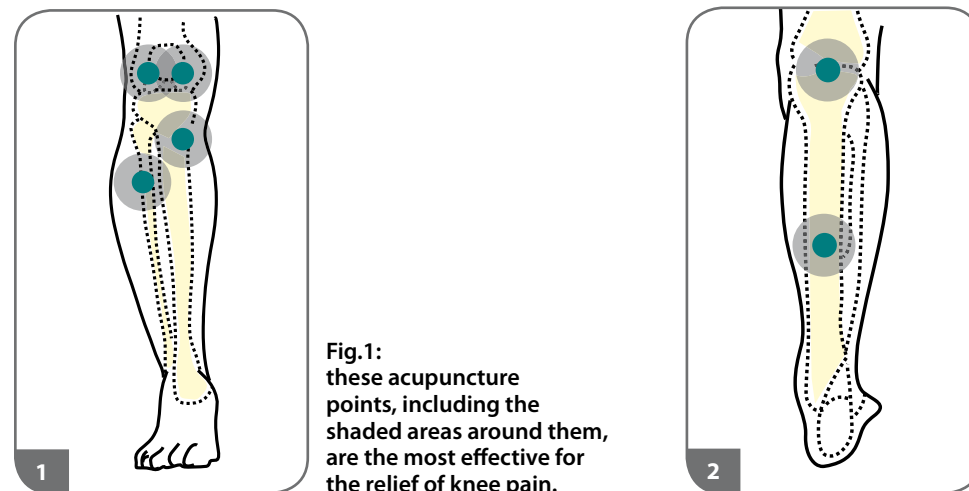
Fig.1:
these acupuncture
points, including the
shaded areas around them,
are the most effective for
the relief of knee pain.

Revitalization - Reenergizing Your Body

People who work late, students who study through the night, and many others suffering from loss of energy of one type or another, will find the revitalizing and invigorating effect of an **ELECTRO-ACUPUNCTURE** application a very valuable asset. Using a circular motion, stimulate the points shown in the drawings with the **ELECTRO-ACUPUNCTURE** device, to achieve a revitalizing and invigorating effect and to enhance the body's energy level, carefully avoiding the areas and conditions noted above in the **WARNING** section.

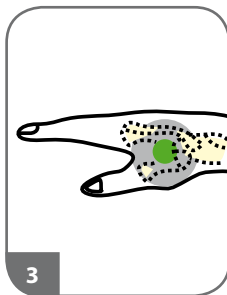
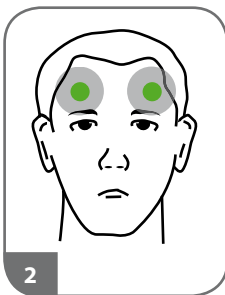
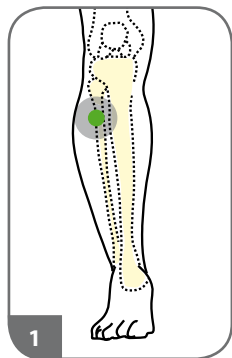
Fig.2 is the most effective in this category.

Fig.3 is also considered important.

In this application, it is very effective to alternate between the different points in fig. 1, 2 and 3.

Classic Acupuncture Points For Specific Areas Of Pain

The user may refer to more complete acupuncture handbooks to add to the application repertoire. In any event, the device can always be applied directly to the general area of the pain, ache or tension, without any consideration of the classic acupuncture point. The **ELECTRO-ACUPUNCTURE** device can only be used for the symptomatic relief of all the previously described conditions, and not as a cure to the cause.

CARE OF YOUR ELECTRO-ACUPUNCTURE DEVICE The **ELECTRO-ACUPUNCTURE** is not only easy to operate, it is simple to maintain. It requires only occasional cleaning with a moist cloth. Do not wet or dip in water. The device operates on a standard 9 volt alkaline battery.


IMPORTANT: The **ELECTRO-ACUPUNCTURE** is only a device to help you in the management of pain and tension. It is by no means a cure. Please follow with attention the instructions, warnings, restrictions and precautions mentioned in this booklet.

BATTERY INSERTION

If battery is inserted correctly the red light on the side of the **ELECTRO-ACUPUNCTURE** will flash when turned on. If battery is inserted incorrectly the red light will not flash but will remain lit when turned on. Turn battery 180 degrees reinsert and check for flashing light.

GUARANTEE CERTIFICATE

Dealer & stamp

Purchase Date

ELECTRO-ACUPUNCTURE is guaranteed for 3 period of one year from the date of original purchase. Should during this period of guarantee the product prove defective due to improper materials or workmanship, the Retailer, Authorised Service Centre or Authorised Servicing Dealer will, without charge for labour or parts, repair or, at its option, replace the product or its defective parts according to the following conditions:

CONDITIONS:

1. This guarantee will be granted only when the original invoice or sales receipt (indicating the date of purchase, product type and dealer's name) is presented together with the defective product.
2. This guarantee will not reimburse nor cover the damages resulting from normal wear or tear, accidents or any misuse of the product nor from any adjustments made to the unit or repair of the product carried out by an unauthorised service centre.
3. This guarantee does not cover transport costs, home service transport costs or any risks of transport relating directly or indirectly to the guarantee of the product.

INFORMATION FOR USERS

In accordance with European Directives 2002/95/CE, 2002/96/CE and 2003/108/CE on the restriction of the use of dangerous substances in electric and electronic equipment as well as their waste disposal. The barred symbol of the rubbish bin shown on the equipment indicates that, at the end of its useful life, the product must be collected separately from other waste.

Do not dispose of electrical appliances as unsorted municipal waste, use separate collection facilities; therefore, any products that have reached the end of their useful life must be given to waste disposal centres specialising in separate collection of waste electrical and electronic equipment, or given back to the retailer at the time of purchasing new similar equipment, on a one for one basis. The adequate separate collection for the subsequent start-up of the equipment sent to be recycled, treated and disposal of in an environmentally compatible way contributes to preventing possible negative effects on the environment and health and optimises the recycling and reuse of components making up the apparatus.

If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being. Contact your local government for information regarding the collection systems available. Abusive disposal of the product by the user involves application of the administrative sanctions according to the laws in force.

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